



CLEANI with Tank Plus

High-Volume Washer for Dental 3D Printing



Ideal for High-Volume Production

Thanks to its high capacity and fast washing cycle, Tank Plus makes CLEANI ideal for demanding digital dentistry workflows, accommodating 8 full arches efficiently.



Made to Match Our New 3D Printing Solution Plus

With Tank Plus, our alcohol washer becomes the perfect fit for our new 3D Printing Solution Plus - wash up to 8 full-arch applications in one go to match the high production that SOL Plus and CURIE Plus deliver perfectly.



Enjoy a Short 5-Min Wash Cycle

Complete an entire washing cycle in just 5 minutes. With Tank Plus, enjoy the advantage of swiftly washing high-volume 3D printed applications in record time.

CLEANI with **Dual Tanks**

Dual Tank Washer for Dental 3D Printing



Perfect for Biocompatible Workflows - Up to 4 Full-Arches at Once

The versatile two-stage cleaning process guarantees a safe and hygienic wash for your 3D-printed applications, such as crowns, bridges, splints, surgical guides, denture bases, teeth, and models, among many others. It is ideal for both biocompatible and non-biocompatible rigorous workflows.



2-in-1: Clean Both Your Build Platform and Applications

Save time and up to 50% alcohol consumption with Dual Tanks. Accommodate a wide range of build platforms with the 3D prints still attached for a simultaneous wash process.



6 Minutes, 2 Cycles

Remove up to 90% of excess resin during the **dirty wash** and employ the second tank for a **clean wash** ensuring the removal of all remaining resin particles.



A Seamless Move from Printer to Washer CLEANI with Dual Tanks & CLEANI with Tank Plus



Perfectly Sized for Dental Practices



No guesswork: Hydrometer for Clean Alcohol in the Tank



3 Pre-set cycles: 3, 5, & 10 mins.



Vortex motion with gentle vibrations to reduce damage risk



23 cm 9.1 in

High

31 cm 12.2 in

24 cm

9.4 in


Width


23 cm 9.1 in

Printing Your Perfect Smile



Technical Specifications

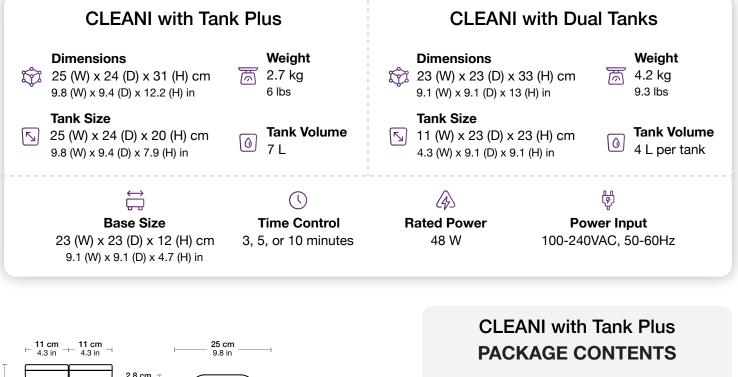

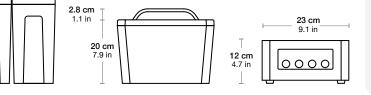

23 cm

9.1 in

High

33 cm 13 in

Depth 23 cm 9.1 in 1 x Base

- 1 x Hydrometer
 - 1 x Power Cable
 - 1 x AC Adapter
- 1 x Mesh Basket

1 x Tank

1 x Lid

CLEANI with Dual Tanks **PACKAGE CONTENTS**

1 x Base	2 x Jig	2 x Hydrometer
2 x Tank	2 x Bar	1 x Power Cable
2 x Lid	2 x Mesh Basket	1 x AC Adapter

Width 23 cm 9.1 in 23 cm 9.1 in Visit ackuretta.com to learn more For sales inquiries, contact sales@ackuretta.com

Depth

